

SRT Model: YC-KIII-132




## Features:

- 1. Carbon fiber material;
- 2. It has one bow form ruler;
- 3. Can be matched with TPS system;
- 4. Comes with one carbon fiber head support and two vacuum fixation mat blockers;
- 5. Compatible with u-shaped head mask, 3-point head mask, 5-point head-neck-shoulders mask, 2-point breast mask, 4-point belly-pelvis mask and 6-point pelvis mask;
- 6. This product can be added metal coordinates into the baseplate, and the central coordinates of the tumor can be read directly through the supporting software, without the need for CTSIM;
- 7. Positional accuracy ≤ 0.5mm.




## Accessories :

Carbon fiber baseplate	1pc	Level governor (screw)	1set
Bridge with laser	1 set	Supine Head Supports	1set
Locating rod slider	2sets	Shell frame for head positioning	1set
Locating rod	2pcs	Vacuum fixation dam-board	1set
Hand grips	2pcs	Vacuum fixation mat repositor	1set

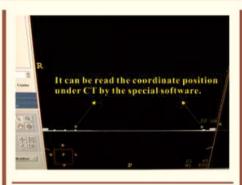




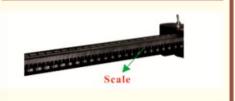








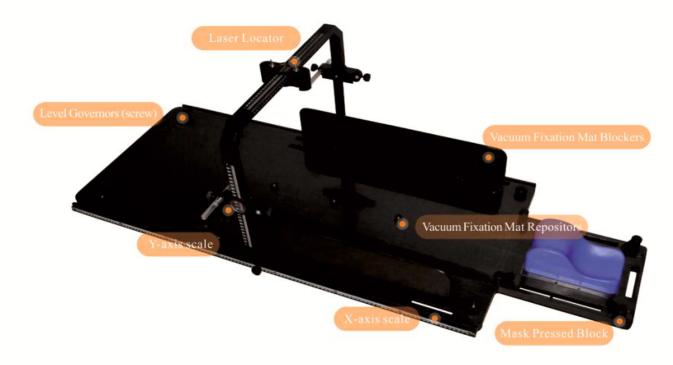

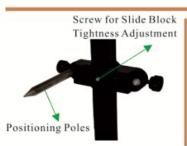








Vacuum Fixation Mat Blockers Model: YC-DF157

- 0 Radiotherapy Positioning Products
- 1. Carbon fiber material;
- 2. It has one bridge with laser lamp;
- 3. Can be matched with TPS system;

4. Maximum patient load: 150kg;

5. Comes with two vacuum fixationmat blockers.

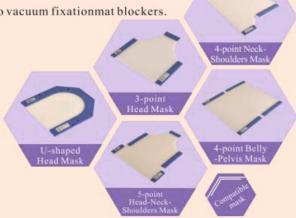






Alloy aluminum material; One set including 2 pieces; This product is mainly used for horizontal adjustment for the baseplate or the CT couchtop.

## Features:

- 1. Carbon fiber material;
- It has one bridge with laser lamp;
- Can be matched with TPS system;
- Comes with two vacuum fixation mat blockers.

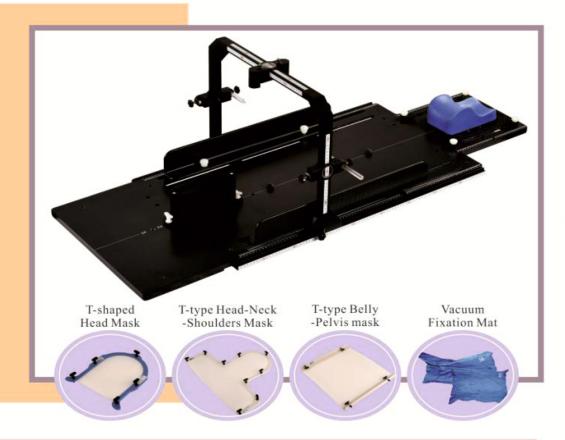

## Vacuum Fixation Mat







- $1. \ Surface \ of the \ product \ is \ comfortable \ and \ light \ with \ good \ gas \ tightness;$
- 2. Used for radiation therapy in patients position and fixation;
- 3. The product filling with EPS grains and the surface is environmental fabric;
- 4. It is automatic sealing, and easy to operate;
- 5. The black protective cap of the product is on the front end, it can prevent the air leakage in the handling process, and it also can have the effect of secondary seal;
- 6. Can customize the card for recording the patient's information;
- 7. Can customize any kinds of sizes;
- 8. Can provide the special vac-lock cushion for Gamma knife.

